

PTHC155GCSSLHD



71.5" x 23.63" x 23.75" (H x W x D)

Single chamber passive heating blanket warmer with a self-closing glass door, stainless steel interior, digital thermostat, stainless steel cabinet, and lock

Highlights:


- Designed and assembled in the USA
- No internal fans or moving parts ensure silent, safe, and efficient operation
- Stainless steel interior ensures durable use and easy cleaning

Product Features:

90 to 140°F temperature range	Designed and set to keep blankets warm in hospitals, medical and physical therapy offices, spas, and other facilities
Passive heating	Fixed heaters located at the bottom of the cabinet with no fans or moving parts to minimize lint buildup for safer, more efficient use
Digital thermostat	Electronic controls are conveniently located for easy temperature adjustment
Electronic readout	Large LED temperature display ensures your staff can easily monitor the interior temperature
Silent operation	No moving parts to ensure a quiet performance
Double pane tempered glass door	Heat-safe glass door offers a full view of stored contents, with a stainless steel door trim to highlight the interior
Professional handle	Stainless steel handle completes the professional look
Self-closing door	Door will automatically close when left ajar, preventing items from accidentally being left exposed to ambient temperatures
Stainless steel interior	Interior walls are constructed from stainless steel to ensure lasting durability and easy cleaning

PTHC155GCSSLHD Specifications:

Overview	
Nominal Height	72.0" (183 cm)
Height of Cabinet	71.5" (182 cm)
Nominal Width	24.0" (61 cm)
Width	23.63" (60 cm)
Nominal Depth	24.0" (61 cm)
Depth	23.75" (60 cm)
Depth with Handle	25.63" (65 cm)
Depth with door at 90°	46.75" (119 cm)
Capacity	12.17 cu.ft. (345 L)
Door	Glass
Cabinet	Stainless Steel
US Electrical Safety	ETL
Canadian Electrical Safety	ETL-C
Amps	5.5
Voltage/Frequency	115 V AC/60 Hz
Weight	210.0 lbs. (95 kg)
Parts & Labor Warranty	1 Year
UPC	761101074078
Warming Cabinet	
Door Swing	LHD
Reversible	No
Adjustable Shelves	Yes
Shelf Type	Chrome
Shelf Qty	6
Thermostat Type	Digital
Interior Height	57.0" (145 cm)
Interior Width	20.0" (51 cm)
Interior Depth	16.5" (42 cm)
Dolly	Optional
Temperature Range	90 to 140°F

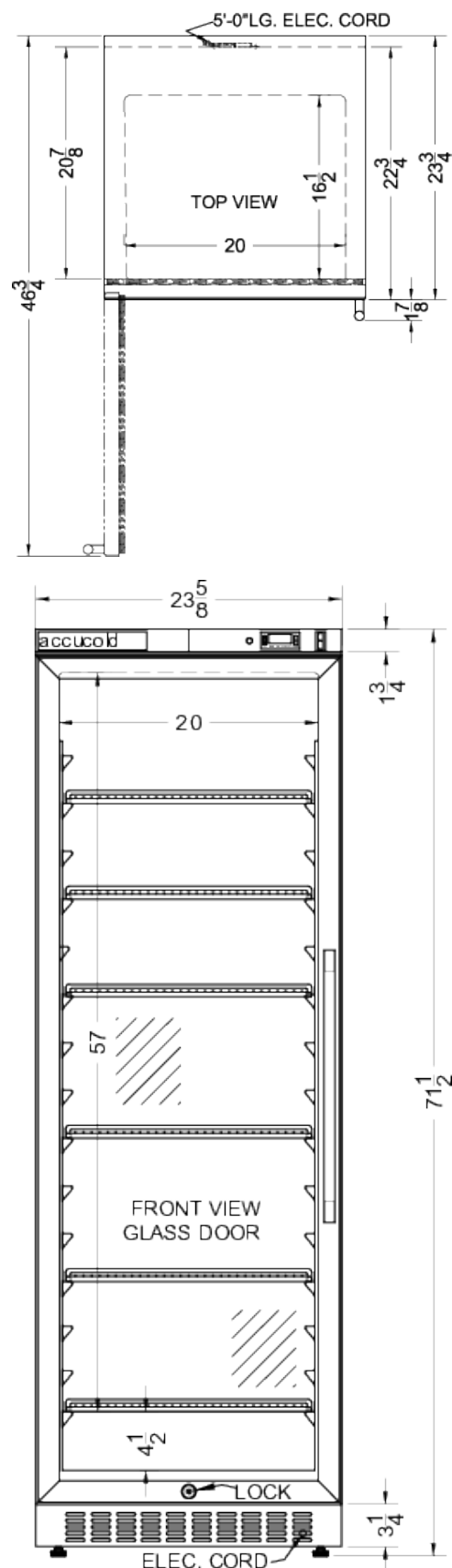

Puretherm Warming Cabinet