

USA Treatment
Tables





Introduction

Chattanooga® is the world's largest manufacturer of rehabilitation equipment for treating musculoskeletal, neurological and soft tissue disorders. For over six decades, Chattanooga has set a benchmark for leadership, reliability and excellence. Our products contribute to better treatment outcomes in hospitals, clinics and home settings worldwide. We lead the physiotherapy industry by example and, through continuous innovation, provide real world solutions for clinicians and their patients.

As the world's leading brand in rehabilitation, Chattanooga is sold in over 80 countries worldwide. Chattanooga is directly represented in our key markets but also has long lasting relationships with dedicated distributors in the USA, Canada, Latin America, EMEA and Asia Pacific. Our Chattanooga distributors and customers are supported for order placement, training, product delivery and after sales service by a dedicated team of highly skilled DJO® employees based all over the world.

Chattanooga is a brand of DJO Global, Inc. DJO Global® provides solutions to aid musculoskeletal health, vascular health and pain management. Our products help prevent injuries or rehabilitate after surgery, injury or from degenerative disease, enabling patients to regain or maintain their natural motion.

Visit www.DJOglobal.com



Contents


Introduction	1
Galaxy® Tables	3
Galaxy® 3 Section	4
Galaxy® 5 Section	5
Montane® Tables	6
Montane® Andes - 7 Section Table	7
Treatment Stools	8
Accessories	9



Treatment Tables

The foundation for any musculoskeletal assessment and treatment is the treatment table. It is also one of the most significant and long standing investments the practitioner will make. The treatment table is one of the foundations of the clinic's style and helps define the atmosphere of the practice. With more than 60 years experience in the rehabilitation arena, Chattanooga® remains the practitioner's partner in design and equipment.

Galaxy[®] Tables

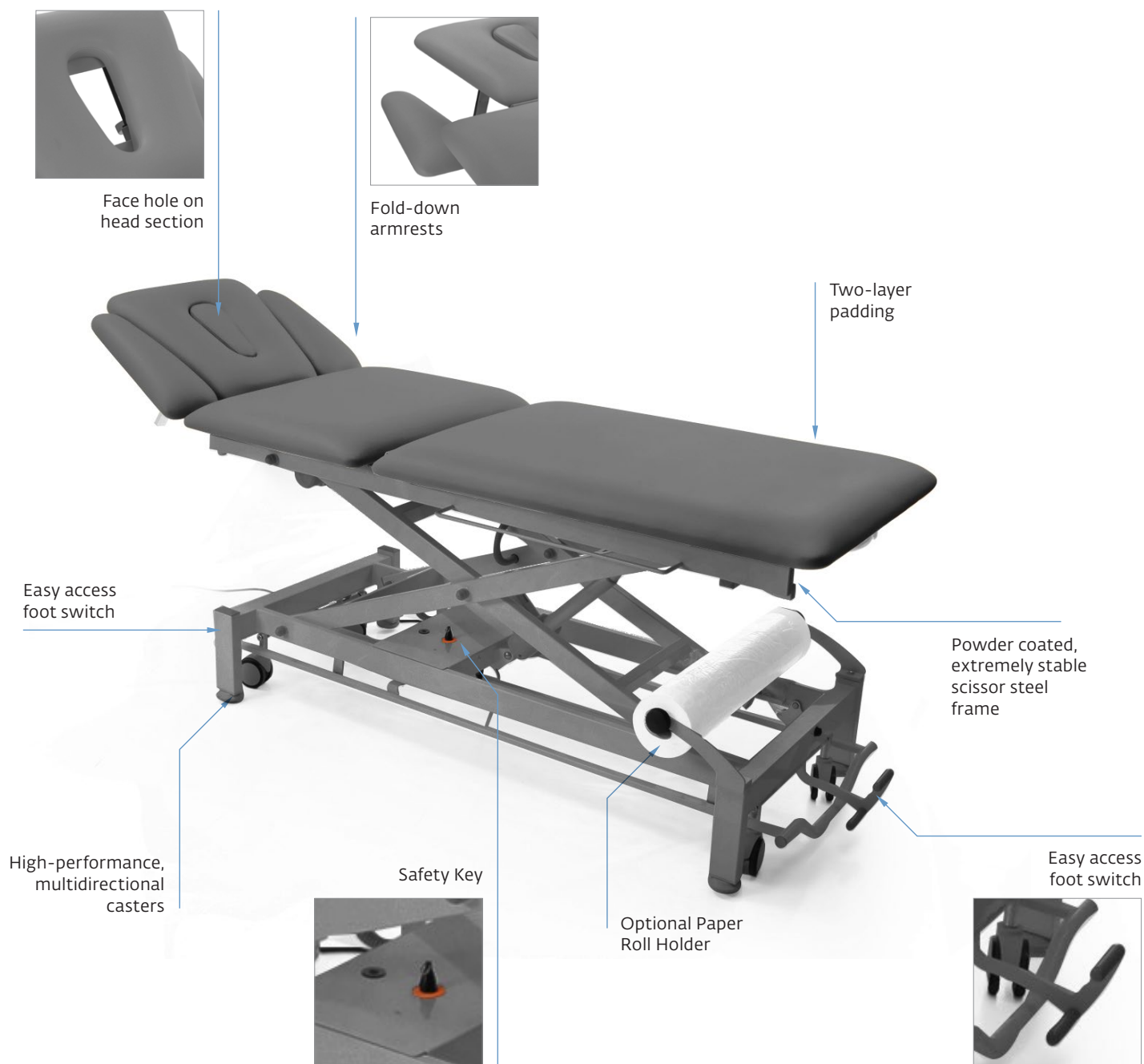

Your table is the core of your practice, and with Chattanooga[®], you have a trusted partner. Every detail of the Galaxy Table range has been designed, engineered and crafted to elegant, exacting standards, delivering a refined experience that maximizes patient comfort while optimizing practitioner access and ease of operation.

Galaxy Tables provide a strong, reliable, stylish and comfortable addition to your clinic. Their slim compact design allows easy access to patients and the scissor lifting frame ensures Galaxy's minimal footprint even when raising patients. Powerful Hallotronic[®] actuators deliver the strength to lift the heaviest of patients, while the optional hydraulic actuation system gives the choice of complete wireless freedom.

Galaxy* 3 Section



Product features:

- Adjustable head section from +40° to -70° with gas spring
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Face hole on head section with removable face hole cover
- Adjustable leg section from 0° to 90° with gas spring
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery
- Powder coated, extremely stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic™ +System electric high/low height adjustment from 21 to 40 inches
- Non-slip rubber feet with 1 cm height adjustment allowing leveling of the table
- 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions Standard



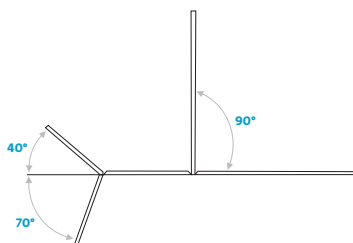


Table Articulations/Dimensions WIDE

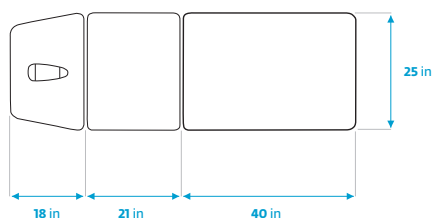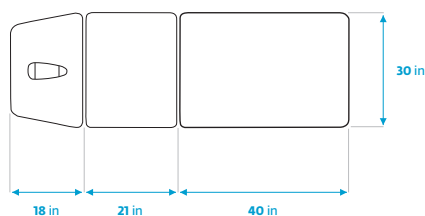

Ordering Information

Part Number Description

Galaxy 3 Section Treatment Table

3133-x	Galaxy 3 Section Electrical, All-round Switch, 4 Castors
3132-x-XL	Galaxy 3 Section Electrical, All-round Switch, 4 Castors, WIDE

Technical Specifications

Electrical

Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Average Weight:	Standard: 274 LBs Wide: 285 LBs
Shipping Weight:	Table Weight + 56 LBs
Table Dimensions*:	Standard: 79 x 25 x 21 - 40 inches Wide: 79 x 30 x 21 - 40 inches
Shipping Dimensions*:	Standard: 86.6 x 33.5 x 31.5 inches Wide: 86.6 x 37.4 x 31.5 inches
Lifting Capacity:	440 lbs Standard 500 lbs Wide
Height Adjustment:	21-40 inches
Electrical safety class:	Class II, Type B

(x) Upholstery Color Options

Beige (x = 101)	Blue (x = 102)	Gray (x = 106)
Graphite Gray (x = 107)	Black (x = 108)	Imperial Blue (x = 119)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.
For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

*= Length x width x height

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO
+Hallotronic™ is a trademark by MEDEN-INMED Spółka z o.o.



Galaxy* 5 Section



Product features:

- Fold-down armrests in head section
- Adjustable head section from +40° to -70° with gas spring
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Face hole on head section with removable face hole cover
- Adjustable leg section from 0° to 90° with gas spring
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery
- Powder coated, extremely stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic™ +System electric high/low height adjustment from 21 to 40 inches
- Non-slip rubber feet with 1 cm height adjustment allowing leveling of the table
- 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

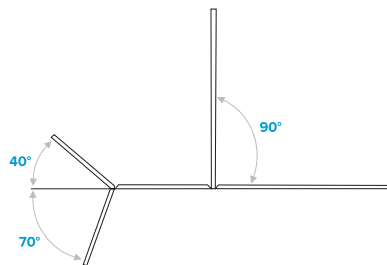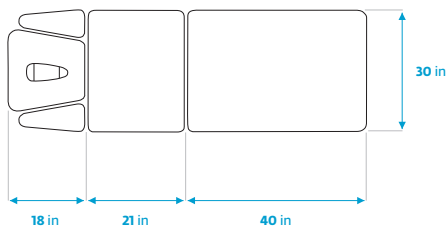

Ordering Information

Part Number Description

Galaxy* 5 Section Treatment Table

3153-x-XL	Galaxy 5 Section Electrical, All-round Switch, 4 Castors, WIDE with Postural Flex
-----------	---

Technical Specifications

Electrical

Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Average Weight:	285 LBS
Shipping Weight:	Table Weight + 56 LBS
Table Dimensions*:	Wide: 79 x 30 x 21 - 40 inches
Shipping Dimensions*:	Wide: 86.6 x 37.4 x 31.5 inches
Lifting Capacity:	500 lbs
Height Adjustment:	21-40 inches
Electrical safety class:	Class II, Type B

Upholstery Color Options

Beige (x = 101)	Blue (x = 102)	Gray (x = 106)
Graphite Gray (x = 107)	Black (x = 108)	Imperial Blue (x = 119)


NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.
For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

* = Length x width x height

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO

+Hallotronic™ is a trademark by MEDEN-INMED Spółka z o.o.

Montane* Tables

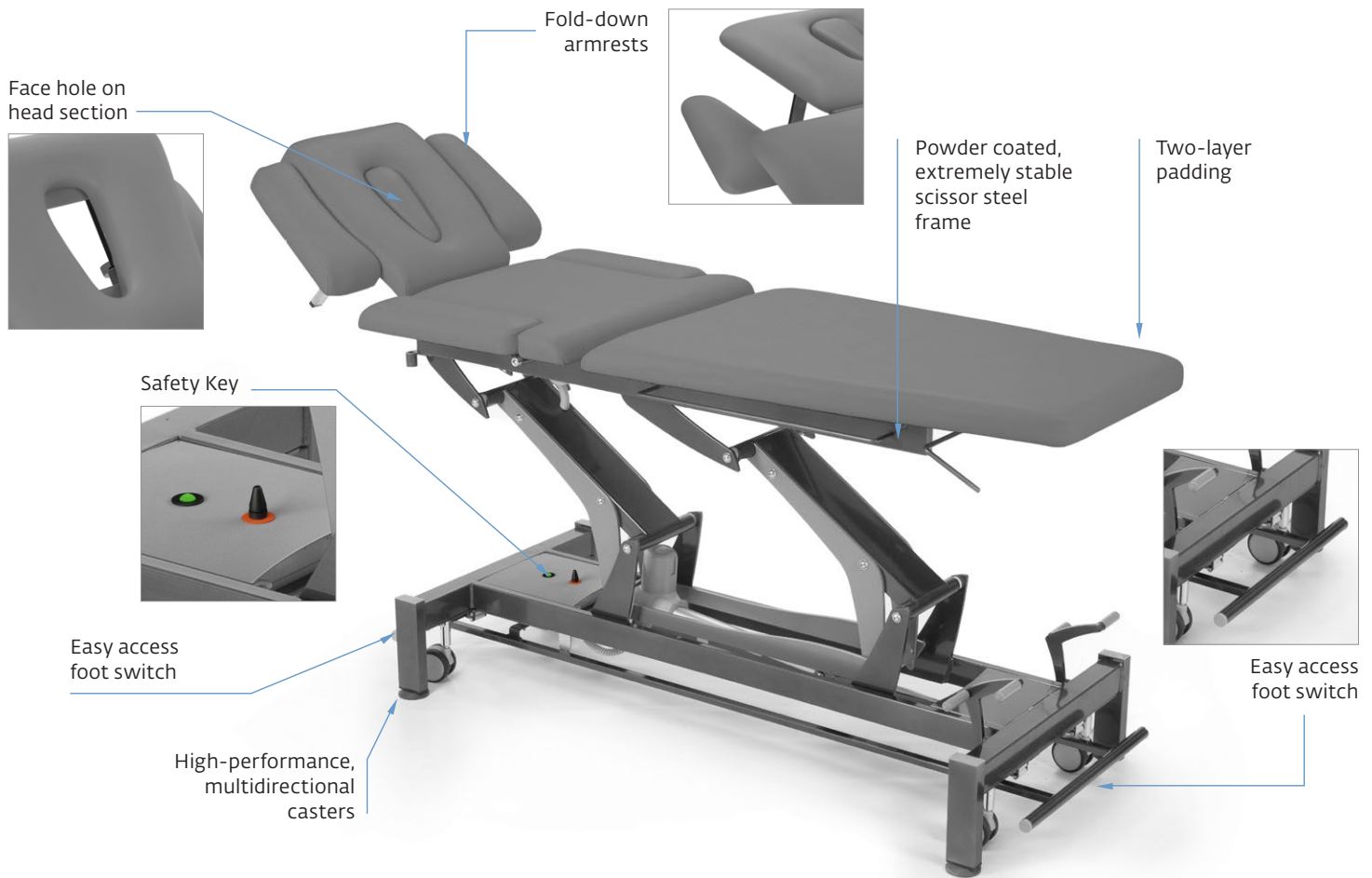

The Montane line is a strong, sturdy, reliable and high quality table that will support your patients and practice for years. Available in a range of color options, Chattanooga® once again offers the practitioner a broad range of real table solutions.

Robust Mechanical Integrity

- Steel frame provides stability during treatment and also during patient transition to and from the table
- Boomerang mechanism provides a smooth vertical lift
- Gas springs assist table adjustment and prevent sudden movements

Quality Design

- Traditional table designed for stability and patient comfort
- A powerful motor provides fast height adjustment
- Easy access foot switch is accessible from any position around the table, making height adjustment easy and convenient
- One movement engages the wheels to move the table or to stabilize it for treatment

Upholstery

Available in 6 fashion colors to complement a variety of settings from contemporary designer trends to classic clinical settings

Easy Table Movement

Includes 4 casters and, with one movement, the wheels engage to move the table or to stabilize it for treatment.

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO

+Hallotronic™ is a trademark by MEDEN-INMED Spółka z o.o.

Montane* Andes - 7 Section Table




Features


- 7 section treatment table includes fold down segments in head and chest section
- Electric high/low height adjustment from 19 to 40 inches for easy and safe patient transfer and adjustment to the ideal working height for treatment
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Features smooth vertical lift with a capacity up to 500 lbs
- Includes 4 casters
- Motorized adjustable middle section with postural flex model
- Adjustable head support - positive/negative inclinations
- Face hole on head section with removable cushion
- Side brackets for stabilizing belts are standard equipment





Table Articulations

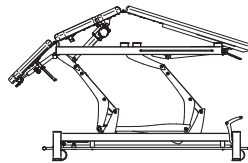

Standard

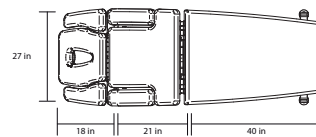

Wide

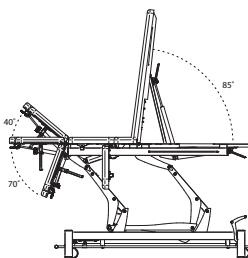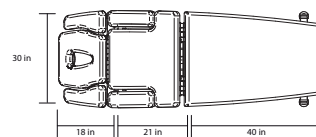

Ordering Information

Part Number Description

Montane Andes

3571-x	Andes 7 Section 120 V - 4 Wheels
3573-x-XL	Andes 7 Section 120 V with Postural Flex - 4 Wheels

Technical Specifications

Height Adjustment:	19 - 40 inches
Mains Power:	230 V / 50 Hz, 120 V / 60 Hz
Standby Consumption:	3.5 W
Input:	120 V / 1.9 A / 60 Hz
Average Weight:	Standard: 274 LBS Wide: 285 LBS
Shipping Weight:	Standard: 330 LBS Wide: 341 LBS
Table Dimensions*:	Standard: 79 x 25 x 21 - 40 inches Wide: 79 x 30 x 21 - 40 inches
Shipping Dimensions*:	Standard: 86.6 x 33.5 x 31.5 inches Wide: 86.6 x 37.4 x 31.5 inches
Lifting Capacity:	440 lbs Standard 500 lbs Wide
Electrical Safety Class:	Class II, Type B
Safety Tests:	Meets Directive 93 / 42 / EEC, Certified to CAN/CSA Standard C22.2 No. 601.1-M90x/A2

* Dimensions are listed as L x W x H

(x) Upholstery Color Options

Beige (x = 01)	Blue (x = 02)	Gray (x = 06)
Graphite Gray (x = 07)	Black (x = 08)	Imperial Blue (x = 19)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO
+Hallotronic™ is a trademark by MEDEN-INMED Spółka z o.o.

Galaxy* / Montane* Stools



Standard Round Stool

Lightweight and comfortable chrome-plated stool for physicians and therapists

- Rubberished castors included
- Round seat diam. 13.7 in
- Height adjustment: 20.5-30 in
- Max. load: 330 LBs
- Base diameter: 19.7 in
- Weight: 13.6 LBs
- Foam density: 60 Kg/m³
- Standard color: Graphite grey
- Upholstery available in all Galaxy/Montane colors



Flexi Saddle Stool


Rubberished castors included

- Tilting seat (front / back): Tilting adjustment 15°
- Height adjustment: 19.7-27.5 in
- Max. load: 330 LBs
- Base diameter: 19.7 in
- Weight: 15.8 LBs
- Foam density: 60 Kg/m³
- Standard color: Graphite grey
- Upholstery available in all Galaxy/Montane colors

Upholstery Color Options						
						
	Beige (x=01)	Blue (x=02)	Gray (x=06)	Graphite Gray (x=07)	Black (x=08)	Imperial Blue (x=19)
STANDARD	371901	371902	371906	371907	371908	371903
FLEXI SADDLE	372001	372002	372006	372007	372008	372003

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO

Accessories

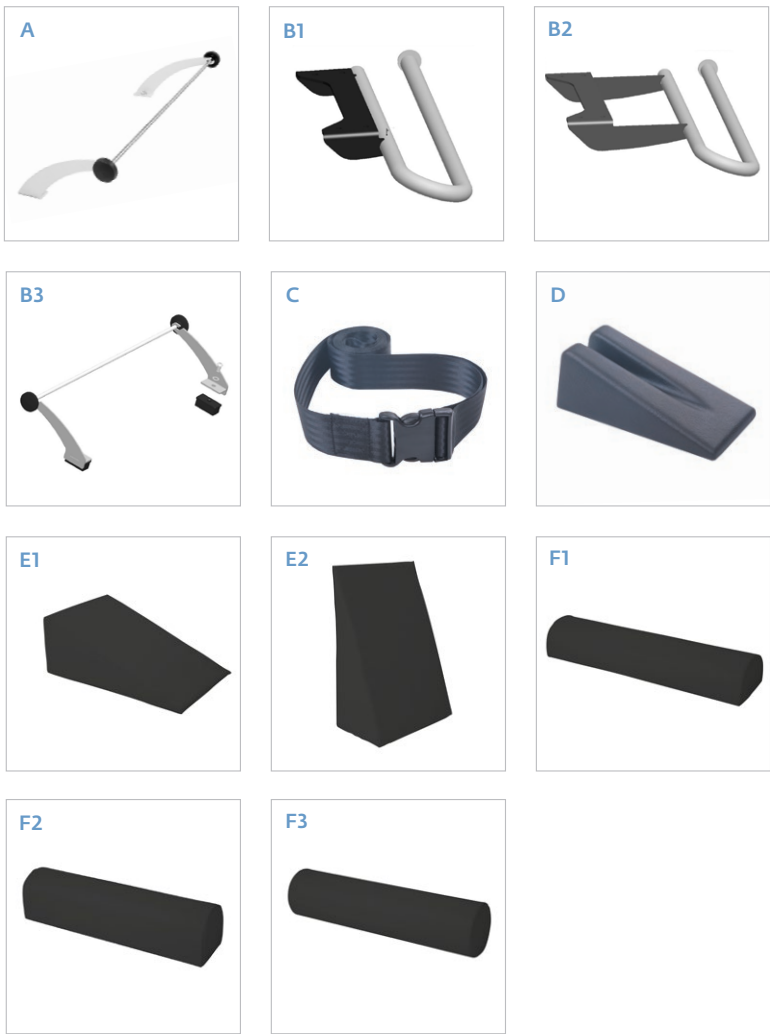

A Paper Roll Holder, Galaxy Tables - Holds paper rolls up to 23.6 inches wide

B Paper Roll Holder, Montane Tables - Holds paper rolls up to 23.6 inches wide

C Stabilization Belt - Available in black, this stabilization belt attaches to side brackets (available on most models)

Wedges and Bolsters are used to reinforce the correct positioning of the patient during treatment.

D Kaltenborn Wedge - For use by professionals trained in manual mobilization techniques.

E1 & **E2** Wedge

F1, **F2** & **F3** Bolster

Ordering Information		
Part Number	Description	
A Paper Roll Holder for Galaxy* Tables		
3110GALXYPPRHLD	Paper Roll Holder for the Galaxy Tables	
B Paper Roll Holders for Montane* Tables		
B1 371008	Standard Montane Tables	
B2 371009	Montane table with Postural Flex	
B3 371010	All Montane tables (universal)	
C Stabilization Belt		
3712	Stabilization Belt 25.8x2 inches (including clasp)	
D Kaltenborn Wedge		
371308	Kaltenborn Wedge 3.7x6x2 inches	
The following items are available in the 6 table colors where 'xx' is replaced by the color of choice.		
E Wedge		
3714-xx	(E1) Wedge 11.8x15.7x4.7 inches	
3715-xx	(E2) Wedge 19.7x13.7x7.9 inches	
F Bolster		
3716-xx	(F1) Bolster 23.6x6x4 inches	
3717-xx	(F2) Bolster 23.6x6x6 inches	
3718-xx	(F3) Bolster 23.6x6x6 inches	
Upholstery Color Options		
Beige (xx = 01)	Blue (xx = 02)	Gray (xx = 06)
Graphite Gray (xx = 07)	Black (xx = 08)	Imperial Blue (xx = 19)
* Dimensions are listed as L x W x H		

* GALAXY is a table name specific to DJO® tables manufactured by MEDEN-INMED exclusively for DJO



DJO, LLC | 1430 Decision Street | Vista | CA 92081-855 | U.S.A.

www.DJOglobal.com

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this catalog do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.

Galaxy and Montane Tables are manufactured by Meden-Inmed Sp. z o.o. exclusively for Chattanooga. A DJO Global Company