


MMT (manual muscle test) kit

User Manual

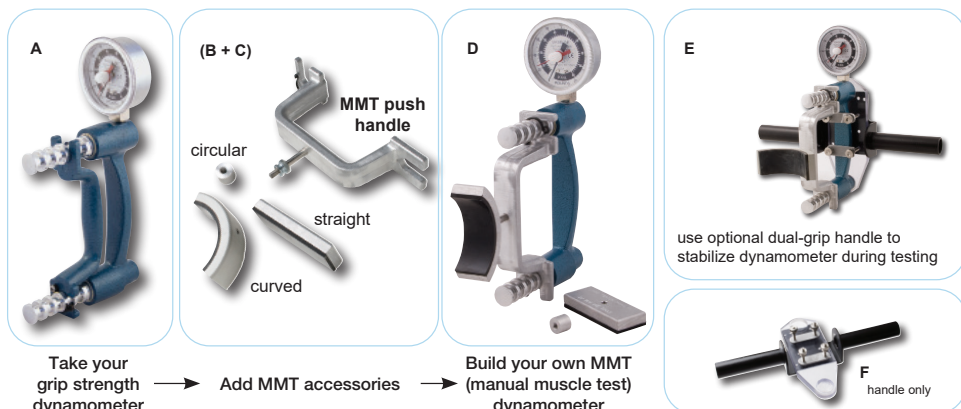

Introduction

Use with any hand dynamometer for MMT testing!

- MMT handle converts ANY hydraulic grip strength (hand) dynamometer to an MMT (Manual Muscle Test) dynamometer
- MMT handle comes with 3 push pad attachments: straight, curved and small circular
- To use simply take a grip strength dynamometer and add your desired MMT accessories to it

WARNING:

- Ensure proper hand position and posture during testing
- Clean grip surfaces between patients
- Do not exceed maximum force capacity
- Stop testing if patient reports pain
- Store in protective case when not in use
- Inspect for damage before each use
- Keep moving parts clean and free of debris
- Check hydraulic system for leaks

Tolerances:

The 200 lb. hydraulic hand dynamometers are tested for accuracy at 20, 60, 100, 140, 180 and 200 lbs. TOLERANCES (Based on 3% through 1st 1/4 of range, 2% through 2nd and 3rd 1/4 of range and 3% through 4th 1/4 of range).

Disposal Method

Dispose of item in accordance with the local/regional/national/international regulations.

