


# single-use iontophoresis electrodes and controllers

## Optima™ electrodes

- pH defining agent adjusts pH to optimize transport when using a carbon conductive element
- pH buffering agent maintains ideal pH
- viscosity-enhancing element creates uniform current distribution for greater patient comfort
- adhesive foam stretches and adheres well when under stress

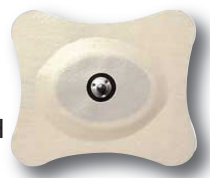

	<b>12 treatment kits per box</b>	
50-0005-0	small 1.5cc	100.00
50-0005-1	medium 2.0cc	100.00
50-0005-2	large 3.0cc	100.00

## Chattanooga® Ionto controller


- dual channel: treat 2 sites or deliver 2 medications simultaneously
- constant current provides precise delivery regardless of skin impedance
- current can be set in 0.1 mA increments between 0.5 mA and 4.0 mA
- automatic 30-second current progression during power on/off for patient comfort
- alert sounds if electrode faults for open current or high impedance, low battery, treatment completion or power left on
- previous dosage can be set as default to maintain frequently-used treatment parameters



00-1335	Ionto controller	385.00
01-3214	lead wire replacements	21.50

## IOGEL electrodes


- clear adhesive allows precise placement on target site
- silver-silver chloride conductor superior to carbon and stabilizes pH without using chemical buffers
- GelSponge® ensures complete contact between the electrode and the skin and provides a conductive medium similar to ultrasound gel
- EZ-Fill Window™ protects adhesive during hydration/handling and no-stick handling

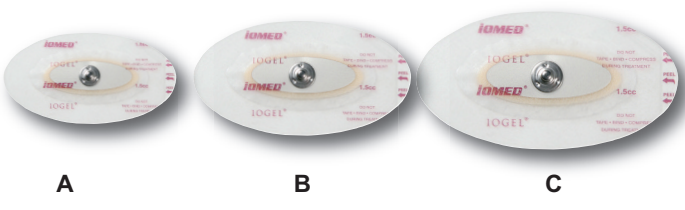

	<b>12 treatment kits per box</b>	
<b>A</b> 50-0002-1	small 1.5cc	100.00
<b>B</b> 50-0002-2	medium 2.5cc	100.00
<b>C</b> 50-0002-3	large 3.5cc	100.00

## Ionto+ Hi-Performance electrodes

- new easy-fill window
- greater adhesive area eliminates pad movement
- silver-silver chloride conductor (superior to carbon) for higher conduction (more drug delivery)
- buffering agents for additional pH stabilization

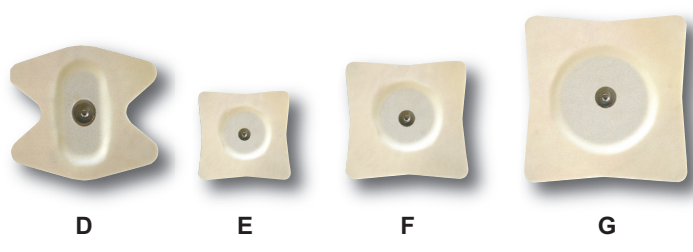

	<b>12 treatment kits per box</b>	
<b>D</b> 13-5250	wings	80.00
<b>E</b> 13-5251	small square	80.00
<b>F</b> 13-5252	medium square	80.00
<b>G</b> 13-5253	large square	80.00

## Trans QE electrodes

- 2 shapes: rectangular and circular for various treatment sites
- adhesive foam stretches and adheres well under stress



	<b>12 treatment kits per box</b>	
50-0000-9	Trans QE small 1.5 - 2.0cc	100.00
50-0000-8	Trans QE medium 2.5 - 3.0cc	100.00

## TransQ Flex electrode

- clover shape: surrounds highly contoured treatment area
- extra adhesive: adds more surface area to increase adhesion when under stress



	<b>12 treatment kits per box</b>	
50-0003-3	TransQ Flex 2.5cc	100.00